



# Ultra Light™ 1 Chin Switch


REVISIONS			
Rev	Description	DATE	BY
1.0	Initial Release	3/6/14	MS, Atec Inc.

## Introduction

The Ultra Light Chin Switch is perfect for activating your favorite device with a simple Chin movement or by opening your mouth.

## Getting Started


- 1) Place the adjustable loop around the neck as shown in the photo to the right.
- 2) Adjust the loop so the end of the switch's lever is where you want it on or near the chin.

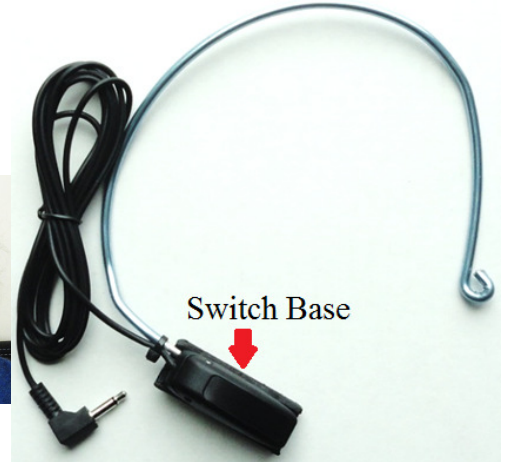

Note: Be sure to take your time and position the switch stably under the chin or the desired location.

## Specifications

Activation Surface Area (in/cm)	1.4" x .5 " (36 x 14 mm)
Activation Force (oz/g)*	< 0.25 oz / < 7.09 g
Travel before Activation (in/cm)*	.09" (2.3mm)
Electrical Rating:	AC 125V 1A
Cable:	6 foot
Connector:	1/8 inch (3.5mm) Right angle mono plug standard for use with PC access, AAC, ECU, adapted Toys, etc.

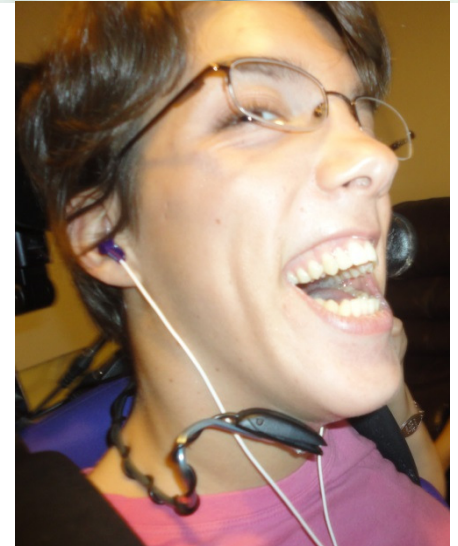

Please don't hesitate to contact us if you have any questions.

Best Regards,  
 Assistive Technology Engineering and Consultation Inc. (Atec Inc.)  
 Sales: +1 (651) 315-7776  
 Fax: +1 (612) 435-2372  
 Email: [Sales@ATengineering.us](mailto:Sales@ATengineering.us)  
 Skype: Sales.atengineering  
 Web: [www.ATengineering.us](http://www.ATengineering.us)